Manual -Laser/torch Mitochondriak® pulse

MANUÁL - Laser Mitochondriak® pulse

MITOCHONDRIAK®

Technické parametre / Technical specifications:

Rozmer laseru / size: 12,5*3*3 cm

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Starting - How it work:

MITOCHONDRIAK®

Who is red light therapy suitable for??

For anyone suffering from body pain, inflammation or needing to improve skin condition.

How to test for photosensitivity?

Before using the torch for red, infrared or blue light therapy, we recommend that you **check your photosensitivity**. To do this, follow these steps: Expose an area, e.g. your forearm, 2-5 cm from the device, to the light for 3 minutes. If you experience any redness or discomfort, you may be photosensitized and we recommend that you consult your doctor or healthcare provider before using the device.

4 phototherapy modes:

1. The first mode is red light therapy (630 + 670 nm).
2. The next mode is near infrared light therapy (NIR = 850 + 900 nm).
3. The third mode is blue light alone, e.g. to reduce inflammation (460 nm).
4. The last mode is all 5 wavelengths shining simultaneously. The mode is a combination of red, NIR and blue light, e.g. to relieve pain and reduce inflammation.

Note: Each of the 4 modes can be activated in 2 modes: the so-called **continuous** mode (rapid pulsation, invisible to the human eye - suitable for regular use and reaching deeper tissues in the body) and in the **10 Hz pulsation mode** (suitable, e.g., for pain relief).

Innovative torch design and features:

The torch features a new pulsation technology that allows the light to penetrate deeper into the skin while preventing burning or overheating the skin.

The torch features an innovative added pulse mode for customers with thin skin to prevent skin burns. It works on the principle of pulsation, which delivers short intermittent pulses to the skin, which makes the energy better absorbed and prevents the skin surface from heating up, while adjusting the rated power.


The red light therapy torch features a 5-minute shut-off timer and is small in size, making it practical and portable. You can use it easily and effortlessly anytime, anywhere on your knees, shoulders, neck, wrists. You can also use it on your fingers, around your waist, ankles or face, or anywhere else you feel pain. This device is also suitable for pets such as dogs, cats, etc..

Thanks to the 10 Hz pulsation mode, the torch can also be a great alternative solution as an analgesic for pain management.

Thanks to its practicality and the fact that it is battery-operated, the torch is an ideal solution for local therapy, with its light penetrating deeper into the body. To reach deeply damaged tissue, the laser actively increases the therapeutic energy using five LEDs with specific wavelengths:

Blue light has a wavelength of 460 nm, **red light** 630 nm and 670 nm, and **near infrared light (NIR)** 850 nm and 900 nm.

Note: Infrared light is invisible to the naked eye. Our laser has 2 infrared diodes (850 and 900 nm), along the edges (see picture).

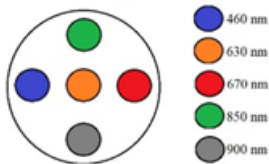How the individual wavelengths work:

Red light penetrates mainly to the skin surface, where it can be absorbed by the mitochondria in the cells.

Improve mitochondrial activity and release more ATP, peroxidase, proteins and other nutrients, thereby increasing cell metabolism and promoting cellular tissue repair.

Red light with a wavelength of 670 nm also improves the viscosity of water in the mitochondria, which will help your mitochondria work better, creating more energy and water.

Red light is therefore suitable for improving the condition of the skin, complexion, collagen, energy and mitochondria..

Infrared light can penetrate into tissues and bones, causing an effect to promote blood circulation, relieve pain and relax muscles. At the same time, thanks to the high intensity and pulsations, the light from the laser penetrates very deeply.

Infrared light is therefore suitable for treating deep parts of the body (bones, joints, tendons, organs, microbiome, ...).

Blue light is now gaining prominence in skin care, including as an **effective acne treatment**. Propionibacterium acnes, the bacteria that causes acne, emits photosensitizers that make the bacteria extremely sensitive to light and susceptible to damage by specific wavelengths. It is blue light that triggers the production of free radicals, which allows these bacteria to be suppressed.

Note: In the Mitochondriak® Pulse torch, we have adjusted the intensity of the blue light, and it only makes up about 2% of the total light. This makes your acne therapy, where you use blue, red and infrared light simultaneously, even more effective and safe.

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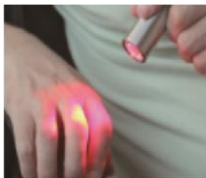
Illumination intensity of mitochondriak® torch pulse:

CONTINUOUS MODE [mW/ cm²]			
Distance	RED	NIR	Všetky
2cm	80-85	90-100	200-220
5cm	40-50	60	110-120
10cm	20	30-35	50-60
20cm	7-8	10-15	20-24

10Hz PULSATION [mW/ cm²]			
Distance	RED	NIR	Všetky
2 cm	40-42,5	45-50	100-110
5 cm	12-15	30	55-60
10 cm	10	15-18	25-30
20 cm	4	5-7,5	10-12

Mitochondriak® pulse torch control:

The laser is controlled by one main button. Each press = another program.

Press 1st time: red (630 + 670 nm) - so-called "continuous mode" is on;

Press 2nd time: red (630 + 670 nm) - pulsation 10 Hz:

Press 3rd time: infrared light (850 + 900 nm) - so-called "continuous mode" is on;

Press 4th time: infrared light (850 + 900 nm) - pulsation 10 Hz:

Press 5th time: blue light (460 nm) - so-called "continuous mode" is on:

Press 6th time: blue light (460 nm) - pulsation 10 Hz;

Press 7th time: all five wavelengths are on - so-called "continuous mode" is on;

Press 8th time: all five wavelengths light up - 10 Hz pulsation:

Press again: laser turns off.

Package contents:

Case

Torch mitochondriak® pulse

Safety glasses

Battery charger

Battery case

Manual + battery

Warnings and Cautions:

Mitochondriak® products, intended for light therapy, are not medical devices. The use of our products is in no way a substitute for a visit to a doctor or medical care.

- Do not look directly into the LEDs or their reflection in the mirror, as this may cause temporary eye irritation.
- We recommend wearing protective glasses (included) when using the device
- Do not use the torch near water.
- Do not use the torch if you are photosensitive.
- Do not allow children to use the torch unsupervised.
- The torch must not be used by sleeping or unconscious persons.
- Do not use the torch on open wounds, as the heat may cause bleeding.
- Do not use the torch in combination with ointments and balms that contain heat-generating ingredients, as this may cause skin burns.
- Using the torch in a manner other than that described in this manual voids the warranty and may result in serious injury.

Warranty conditions:

We guarantee that this equipment will operate within its specifications for a period of 2 years from the date of delivery.

What does this warranty cover??

This limited warranty covers all defects in material or workmanship under normal use as described above during the warranty period. During the warranty period, our company will repair or replace, free of charge, the product or parts of the products that prove to be defective due to faulty material or workmanship under normal use and maintenance.

What is not covered by the warranty?

The limited warranty does not cover any problem that is caused by damage or the result of improper maintenance or self-modification. It also does not cover damage to the product due to natural disasters such as flooding, falling of the product from a height or other natural disasters, theft or loss of the product..

hello@mitochondriak.com

www.mitochondriak.com

instagram/mitochondriak.europe

EasyLight
MITOCHONDRIAK®